# **LIFE-BASE 1 NG**

# **LIFE-BASE 1 NG XL**

# **LIFE-BASE 1 NG XS**

**Portable system / Sistema de transporte**

*Device description and instructions for use*

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
# Contents

---

<b>1. Overview</b>	<b>4</b>	<b>7. Function check</b>	<b>26</b>
1.1 LIFE-BASE 1 NG/ LIFE-BASE 1 NG XL	4	7.1 Visual inspection	26
1.2 LIFE-BASE 1 NG XS	7	7.2 Check catch	26
1.3 Wall mounting	9	7.3 Check protective bag (if present)	26
1.4 Special markings on the device	10	<b>8. Servicing</b>	<b>27</b>
1.5 Safety information in these instructions	11	8.1 Disposal	27
<b>2. Description of device</b>	<b>12</b>	<b>9. Scope of supply</b>	<b>28</b>
2.1 Intended use	12	9.1 General information	28
2.2 Description of function	12	9.2 Standard scope of supply	28
<b>3. Safety instructions</b>	<b>14</b>	9.3 Accessories and spare parts	29
3.1 Safety information	14	<b>10. Technical data</b>	<b>30</b>
<b>4. Assembly</b>	<b>16</b>	<b>11. Warranty</b>	<b>31</b>
4.1 Wall mounting	16	<b>12. Declaration of conformity</b>	<b>32</b>
4.2 Device (not mounted in the case of portable systems with an adapter plate)	16		
4.3 Protective bag (LIFE-BASE 1 NG/ LIFE-BASE 1 NG XL)	16		
4.4 Protective bag (LIFE-BASE 1 NG XS)	17		
4.5 Accessories bag (LIFE-BASE 1 NG/ LIFE-BASE 1 NG XL)	19		
4.6 LIFE-BASE 1 NG XS carrying strap	20		
<b>5. Operation</b>	<b>21</b>		
5.1 Hang portable system in wall mounting	21		
5.2 Hang portable system from stretcher or hospital bed (LIFE- BASE 1 NG/LIFE-BASE 1 NG XL)	22		
5.3 Hang portable system from stretcher or hospital bed (LIFE-BASE 1 NG XS)	23		
5.4 After use	23		
<b>6. Hygiene treatment</b>	<b>24</b>		
6.1 Portable system	24		
6.2 Protective bag and accessories bag	24		
6.3 Attached parts	25		

# 1. Overview

## 1.1 LIFE-BASE 1 NG/LIFE-BASE 1 NG XL
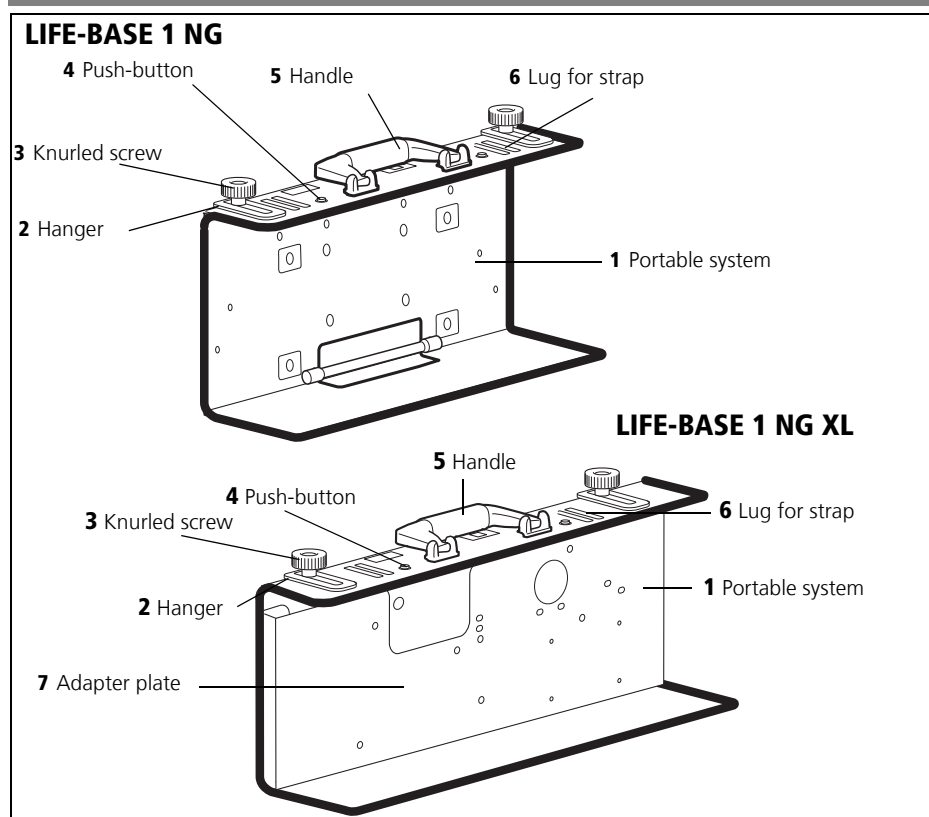
### Key

#### 1 Portable system

The portable system is for storing a combination of ventilation and therapy devices and the necessary components.

#### 2 Hanger

The hanger is for hanging the portable system from a hospital bed or stretcher.

#### 3 Knurled screw

The knurled screw is for fixing the hanger in position.

#### **4 Push-button**

The push-buttons are used to fix the protective bag to the portable system.

#### **5 Handle**

Use the handle to lift and carry the portable system.


#### **6 Lug for strap**

The strap for the portable system can be attached to this lug.


#### **7 Adapter plate**

The devices can be attached to the portable system with the aid of the adapter plate.


**Rear view LIFE-BASE 1 NG**

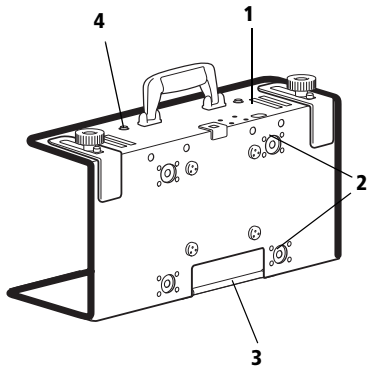**Protective bag**

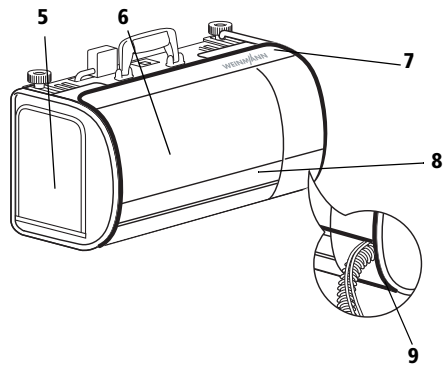**Rear view LIFE-BASE 1 NG XL**

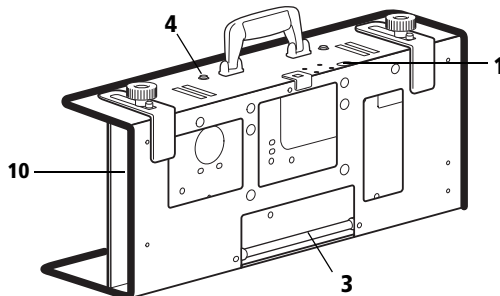
## Key


- |  |  |
|--|--|
| <b>1</b> Cable duct for charger interface                | <b>5</b> Side access compartment                           |
| <b>2</b> Vibration dampers (only with MEDUMAT Transport) | <b>6</b> Viewing window                                    |
| <b>3</b> Hinge for bottom rail                           | <b>7</b> Velcro tab to access the control panel            |
| <b>4</b> Push-button to fix protective bag in position   | <b>8</b> Reflective strip                                  |
|  | <b>9</b> Opening for passing through a patient tube system |
|  | <b>10</b> Adapter plate                                    |

## 1.2 LIFE-BASE 1 NG XS


### Front view

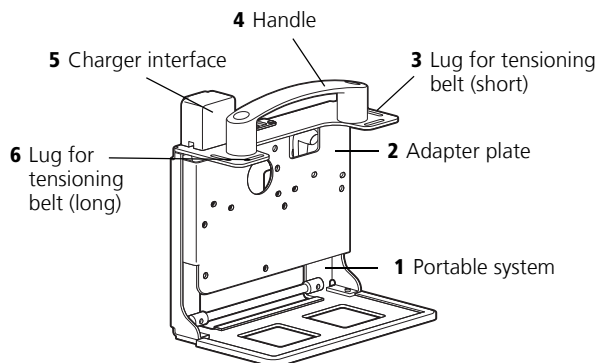
### Rear view

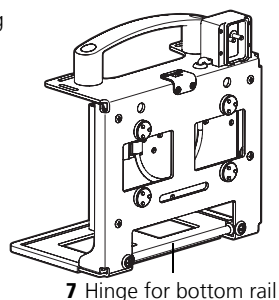
### Protective bag for defibrillator

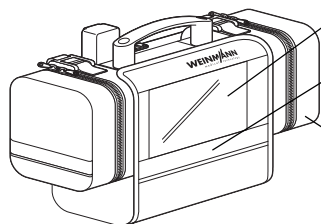
### Protective bag for ventilator
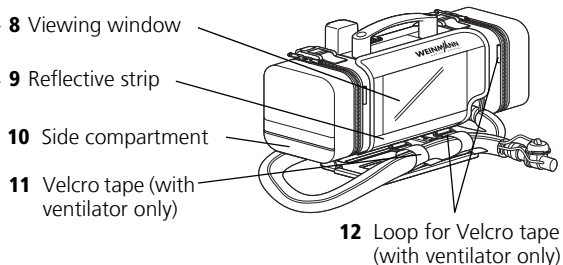
### Carrying strap
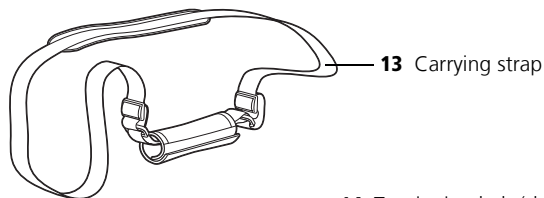
### Tensioning belts
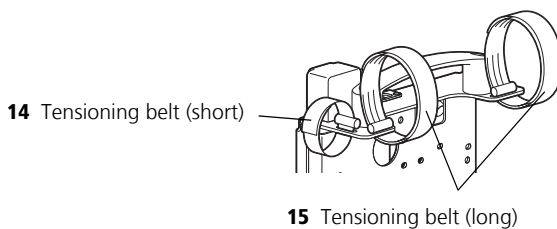
## **1 Portable system**

The portable system is for storing a device and its components.

## **2 Adapter plate**

You can use the adapter plate to attach the device to the portable system.

## **3 Lug for tensioning belt (short)**

You can attach a tensioning belt (short) to this lug.

## **4 Handle**

The portable system can be lifted and carried by the handle. You can attach the carrying strap to the handle.

## **5 Charger interface**

You can charge the device on the portable system using the charger interface.

## **6 Lug for tensioning belt (long)**

You can attach a tensioning belt (long) to this lug.

## **7 Hinge for bottom rail**

This hinge is used to hang the portable system in the wall mounting.

## **8 Viewing window**

This is where to read off the device display when the protective bag is closed.

## **9 Reflective strip**

The reflective strip means the protective bag is also visible in the dark.

## **10 Side compartment**

You can store the components and accessories of the device in the side compartments. With the ventilator, the left-hand side compartment provides access to the compressed gas connection and space for the pressure tube, and the right-hand side compartment provides access to the device inlet filter. With the defibrillator, the left-hand side compartment provides access to the rechargeable battery, and the right-hand side compartment provides space for components such as cables, for example.

## **11 Velcro tape (with ventilator only)**

The Velcro tapes are used to attach the patient tube system to the protective bag.

## **12 Loop for Velcro tape (with ventilator only)**

The Velcro tapes can be attached to these loops which hold the patient tube system.

## **13 Carrying strap**

You can use the carrying strap to carry the portable system on your shoulder.

## **14 Tensioning belt (short)**

The tensioning belt holds the pressure tube of the ventilator if you are not using a protective bag.


## **15 Tensioning belt (long)**

The tensioning belt holds the patient tube system of the ventilator if you are not using a protective bag.


## 1.3 Wall mounting


### BASE-STATION 1 NG
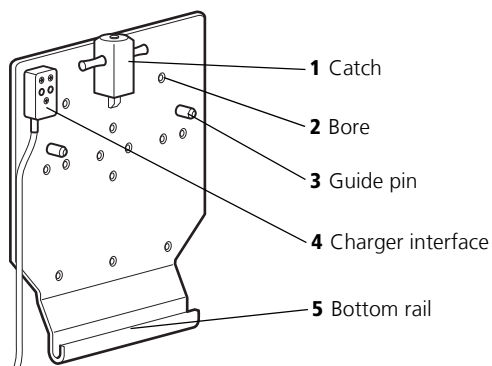
### Catch for wall mounting

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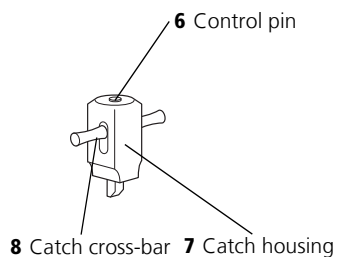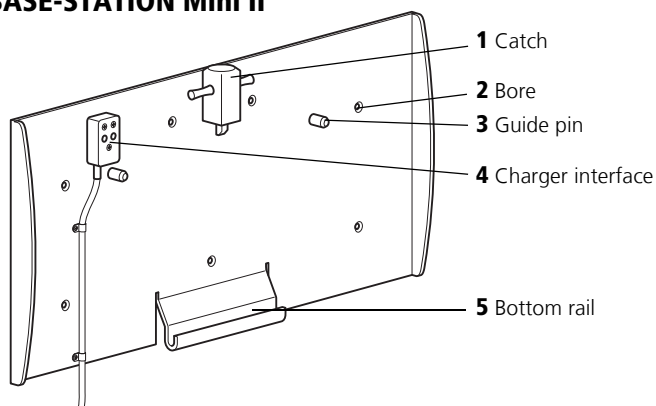**1 Catch**

**2 Bore**

**3 Guide pin**

**4 Charger interface**


**5 Bottom rail**

**6 Control pin**

**7 Catch housing**

**8 Catch cross-bar**

# 1.4 Special markings on the device
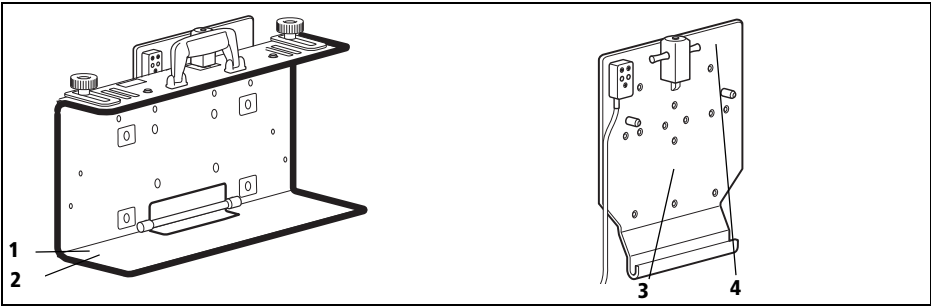
## Markings on the portable system

No.	Symbol	Significance
<b>Information plate for portable system</b>		
1		Year of manufacture
	<b>SN</b>	Serial number of portable system
		CE symbol (confirms that the product conforms to applicable European directives)
<b>Portable system</b>		
2		Do not dispose of the portable system in domestic waste.

## Markings on the wall mounting

No.	Symbol	Significance
<b>Information plate for wall mounting</b>		
3		Year of manufacture
	<b>SN</b>	Serial number of wall mounting
		Do not dispose of the wall mounting with domestic waste.
<b>Wall mounting</b>		
4		Shows how the catch unlatches and the portable system is removed from the wall mounting.

## Markings on the carrying strap

Symbol	Significance
	Do not iron
	Wash at 30 °C
	Do not dry in a tumble-dryer
	Follow instructions for use
	Maximum load

## Markings on the packaging

Symbol	Significance
	Permitted temperature for storage: -40 °C to +70 °C
<b>rh % 0-95</b>	Permitted humidity for storage: max. 95 % relative humidity

## 1.5 Safety information in these instructions

Safety information in these instructions for use is marked as follows:

### **Warning!**

Warns of risk of injury and potential material damage.

### **Caution!**

Warns of material damage and potentially false therapy results.

### **Note**

Contains useful tips.

## 2. Description of device

---

### 2.1 Intended use

---

The portable system is for storing a combination of devices from the MEDUMAT, MODUL and MEDUCORE series, together with the required components (charger interface, tubes, accessories bag, protective bag and electrical cables). This provides a user with mobile units for treating respiratory disorders, for monitoring cardiovascular parameters and for resuscitation support.

The portable system can be latched in the associated BASE-STATION 1 NG and BASE-STATION Mini II wall mountings to fix it securely to the wall.

The BASE-STATION 1 NG can be attached to various retaining devices in emergency vehicles and in hospitals using a variety of attachment sets.

### 2.2 Description of function

---

All the devices provided are securely mounted on the portable system, making them easier for you to transport and operate.

Using the hangers, you can hang the portable system from a stretcher, a hospital bed or a standard hospital rail.

The portable system can be secured in a hospital environment or in an emergency vehicle in conjunction with the BASE-STATION 1 NG or BASE-STATION Mini II wall mounting (see "9.3 Accessories and spare parts" on page 29).

The portable system can also be secured to the ceiling in conjunction with the BASE-STATION 1 NG roof mounting.

The BASE-STATION 1 NG and BASE-STATION Mini II wall mountings are assembled for fixing to the wall. You can also expand them by adding the following attachment sets:

- attachment set for standard hospital rail (BASE-STATION 1 NG and BASE-STATION Mini II):  
attach the wall mounting to a standard hospital rail in vehicles or in hospital.
- attachment set for fitting to a rail (BASE-STATION 1 NG):  
attach the wall mounting to a rail (Ø 19 mm - 38 mm) in hospital. The rail can be horizontally or vertically aligned.

- supplementary attachment set for standard hospital rail (BASE-STATION Mini II):  
for attaching the wall mounting to a standard hospital rail in vehicles

You can remove the portable system from the wall mounting in just one maneuver.

## 3. Safety instructions

---

### 3.1 Safety information

---

Read these instructions for use carefully. It is a component part of the following portable systems:

- LIFE-BASE 1 NG
- LIFE-BASE 1 NG XL
- LIFE-BASE 1 NG XS

The instructions for use must be available at all times. Use the portable system exclusively for the intended purpose described (see “2.1 Intended use” on page 12).

For your own safety and the safety of your patients and in accordance with the requirements of Directive 93/42/EEC, please note the following.

#### Use of the portable system

##### **Warning!**

Use the portable system only if you are medically trained and have received instruction in ventilation/aspiration. Severe physical damage may be caused by incorrect use.

##### **Caution!**

- Before working with the portable system, you must have understood how to use it.
- To avoid infection or bacterial contamination, follow the section entitled “6. Hygiene treatment” on page 24.
- Do not use the portable system in the vicinity of a magnetic resonance tomography (MRT) device.
- Use the portable system only for the intended use described (see “2.1 Intended use” on page 12).
- During the journey/flight, the portable system must be hung in the wall mounting intended for it, as otherwise it might fall over or be thrown around the vehicle. The wall mounting must be securely attached. Follow the fitting instructions for the wall mounting and the attachment set in question.
- Note that the portable system may not be loaded with an additional weight, e.g.

be leant on by a person, even once it is fixed in position. The portable system may fall down and injure people or the device may be damaged.

## Use of the wall mounting

### Warning!

- Only the portable systems LIFE-BASE 1 NG, LIFE-BASE 1 NG XL, LIFE-BASE 1 NG XS and LIFE-BASE Mini II should be suspended in the BASE-STATION 1 NG and BASE-STATION Mini II wall mountings. The weight and dimensions of these portable systems are exactly matched to the wall mountings.
- Hanging third-party products or adapters from other manufacturers is not permitted, as this does not comply with the use in accordance with purpose of this wall mounting.

## Accessories/replacement parts

### Caution!

If third-party items are used, functional failures may occur and fitness for use may be restricted. Biocompatibility requirements may also not be met. Please note that in such cases, any claim under warranty and liability will be voided if neither the accessories nor genuine replacement parts recommended in the instructions for use are used.

## Maintenance

### Caution!

Have servicing and maintenance work carried out only by the manufacturer, WEINMANN Emergency, or by trained staff.

## 4. Assembly

---

### 4.1 Wall mounting

---

The wall mounting for the portable system is ready for fitting to a wall. Using a suitable attachment set (see Section 9.3 on page 29), you can also fit the wall mounting to a standard hospital rail or an ordinary rail (Ø 19 mm - 38 mm). When fitting, follow the relevant fitting instructions enclosed with the wall mounting and the attachment sets in question.

### 4.2 Device (not mounted in the case of portable systems with an adapter plate)

---

**Note:**


In most cases, the device is already mounted on the portable system. You only need to mount the device if you order a device separately.

1. Unscrew the adapter plate from the portable system.
2. Screw the adapter plate onto the rear of the device.
3. Screw the adapter plate and the device to the portable system.

### 4.3 Protective bag (LIFE-BASE 1 NG/ LIFE-BASE 1 NG XL)


---

#### Attach protective bag

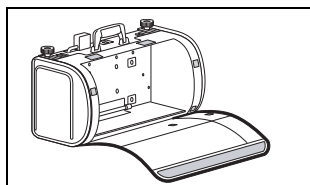
1. Place the protective bag on a firm, level surface.
2. Place the portable system inside, on the base of the protective bag.
3. Attach the protective bag to the portable system using the push-buttons.

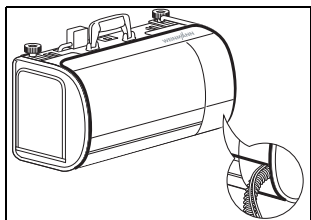
4. If there is a ventilator on the portable system: Route the patient tube system through the front opening in the protective bag.
5. Close the velcro tabs on the front and rear.

## Remove protective bag

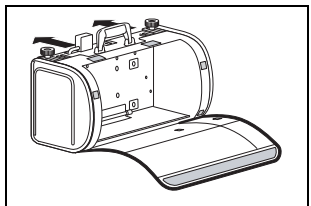
1. Place the portable system on a firm, level surface.
2. Open the Velcro tabs on the front.
3. Open the push-buttons on the portable system and remove the portable system from the protective bag from the rear.

## 4.4 Protective bag (LIFE-BASE 1 NG XS)


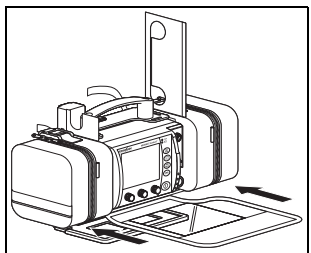
1. Place the portable system on a firm, level surface.
2. Remove all the device components.
3. On the ventilator: push the underside of the protective bag into the space between the portable system and the ventilator from the front and pull the top of the protective bag over the portable system as far as the handle.

### Alternatively

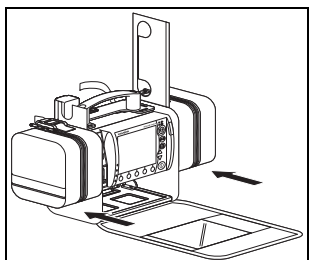
on the defibrillator: push the protective bag over the portable system from the front.

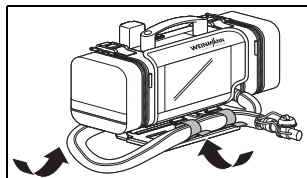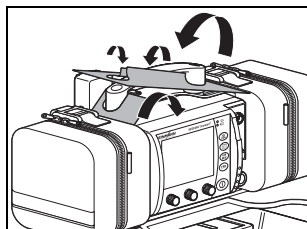
4. Close the Velcro tab on the rear of the portable system.
5. Place the protective bag around the handle and close the Velcro tabs under and next to the handle.
6. Connect all components to the device.
7. On the defibrillator: stow the electrodes under the defibrillator.
8. Close the viewing window of the protective bag.
9. On the ventilator: clamp the patient tube system in the space between the portable system and the ventilator and fix it in position with the Velcro tabs on the protective bag.
10. Stow the accessories and small parts in the side compartments of the protective bag.

## Remove protective bag

1. Place the portable system on a firm, level surface.
2. On the ventilator: open the Velcro tapes on the patient tube system and pull it out of the space between the portable system and the ventilator.
3. On the ventilator: open the viewing window of the protective bag.
4. On the ventilator: detach the patient tube system.
5. Open the Velcro tab on the rear of the portable system.

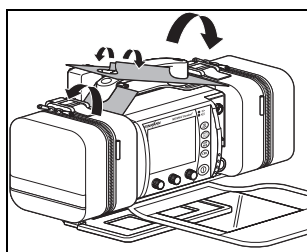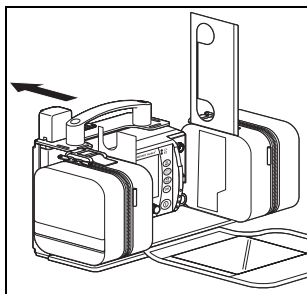6. Take the portable system out of the protective bag backwards.


## 4.5 Accessories bag (LIFE-BASE 1 NG/ LIFE-BASE 1 NG XL)

---

There are four types of accessories bag:


- accessories bag for ventilators (fitted on the side)
- accessories bag for defibrillation/ventilation (fitted on the side)
- accessories bag for defibrillation (fitted on the portable system)

### Fit accessories bag for ventilators

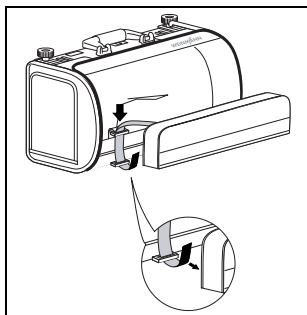1. Open the Velcro tabs on the rear of the accessories bag.
2. Pass the Velcro tabs through the loops on the protective bag.
3. Close the Velcro tabs on the rear of the accessories bag.

## Fit accessories bag for defibrillation/ventilation

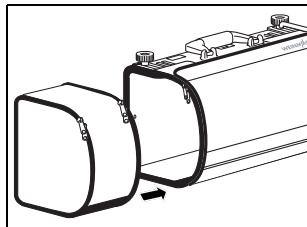
1. Hook the zip on the accessories bag into the zip on the protective bag.
2. Close the zip.

### **Note**

The accessories bag for defibrillation is attached to the left of the protective bag, the accessories bag for ventilation to the right of the protective bag.


## Fit accessories bag for defibrillation

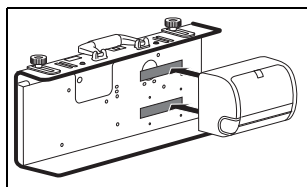1. Press the accessories bag against the rear wall of the portable system.

## 4.6 LIFE-BASE 1 NG XS carrying strap

---

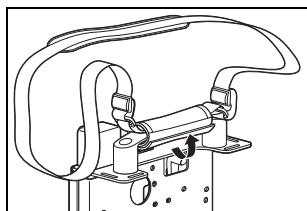1. Place the Velcro tab of the carrying strap around the handle and close it.

# 5. Operation

---

## 5.1 Hang portable system in wall mounting


---

### Hang portable system

#### Note

If you are using a protective bag, ensure that the protective bag is properly located first.

1. To hang the portable system in the wall mounting, place the portable system in the bottom rail.


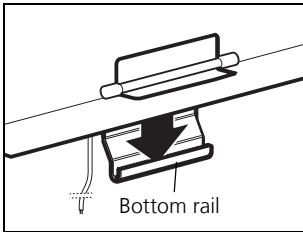
#### Warning!


#### Risk of injury from loose portable system!

An inadequately secured portable system may fall off or be thrown around the vehicle.

- Ensure that you both hear and feel the catch engage.
  - Check that the **red control pin** is flush with the top edge of the catch housing and does not protrude.
2. Swivel the portable system towards the wall until you hear and feel the catch cross-bar engage.
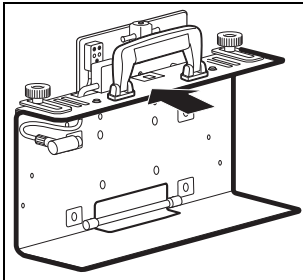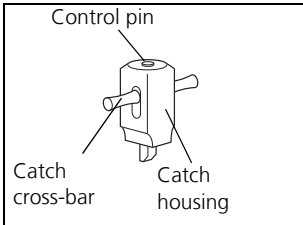
## Remove portable system

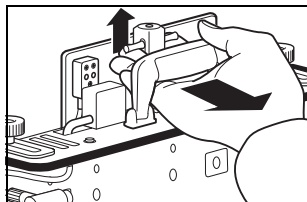1. For unlatching, it is essential to grasp **through** the handle so that the portable system does not fall out of the wall mounting and injure someone.
2. To unlatch, push the catch cross-bar upwards.
3. Swivel the portable system forwards a little way and lift the portable system up off the bottom rail.

## 5.2 Hang portable system from stretcher or hospital bed (LIFE-BASE 1 NG/ LIFE-BASE 1 NG XL)


---

You can also hang the portable system from a stretcher or a hospital bed, as long as the rail diameter does not exceed 35 mm.

### Caution!

Note that the portable system may not be loaded with an additional weight, e.g. be leant on by a person, even once it is fixed in position. The portable system may fall down and injure people or the device may be damaged.

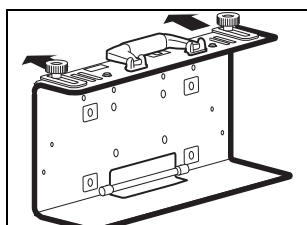1. Undo the knurled screws.
2. Extend the hanger the required distance.
3. Tighten the knurled screws back up.
4. Hang the portable system from a stretcher, a hospital bed or a standard rail.

### Note

The portable system can only be hung back on a wall mounting once you have pushed the hanger back into its original position.

## 5.3 Hang portable system from stretcher or hospital bed (LIFE-BASE 1 NG XS)


---

You can only hang the portable system from a stretcher or a hospital bed if you use a protective bag.

### Caution!

Note that the portable system may not be loaded with an additional weight, e.g. be leant on by a person, even once it is fixed in position. The portable system may fall down and injure people or the device may be damaged.

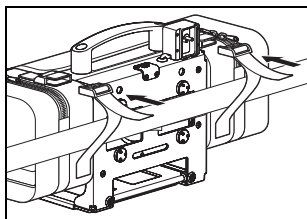1. Open the magnetic closures on the straps on the rear of the protective bag sideways.
2. Guide the straps round a rail of the stretcher or hospital bed.
3. Close the magnetic closures.

## 5.4 After use

---

1. Open the protective bag (if present).
2. If present: Remove the pressure tube of the oxygen supply or stow it in the side compartment of the protective bag.
3. If present: Remove the patient tube system.
4. Close the protective bag again (if present).

# 6. Hygiene treatment

The portable system must be cleaned and disinfected after every use.  
Always perform a function check afterwards (see "7. Function check" on page 26).

## 6.1 Portable system

The portable system is cleaned without its devices by the scouring and wiping method. You can take all the devices off the portable system for this purpose.

- See the instructions for use for the disinfectant used.
- For disinfecting by wiping, we recommend terralin® protect.
- We recommend wearing suitable gloves (e.g. household or disposable gloves) for disinfecting.

**Note**  
For further information on hygiene treatment and a list of all cleaning agents and disinfectants which can be used, please see our Internet brochure at [www.weinmann-emergency.de](http://www.weinmann-emergency.de).

## 6.2 Protective bag and accessories bag

Subject the bags to hygiene treatment in accordance with the table below:

Bag	Cleaning	Disinfecting	Sterilizing
PVC protective bag PVC accessories bag	In hot water with a mild detergent	Disinfectant wipe	Not permitted
CORDURA® protective bag CORDURA® accessories bag	In hot water with a mild detergent <b>Alternatively:</b> 30 °C cycle in a washing machine, no spin	30 °C cycle in a washing machine (no spin), using a suitable disinfecting additive (disinfecting detergent or impregnation agent)	Not permitted

**Note:** To prevent unattractive washing powder deposits, we recommend using a liquid detergent.

## 6.3 Attached parts

---

Follow the instructions for cleaning, disinfecting and, if appropriate, sterilizing in the associated instructions for use.

# 7. Function check

---

You must subject the portable system to a function check before it is used each time and every time it has been dismantled, but in any case at least every 6 months. If you discover faults during the function check, you may not use the portable system again until the faults have been rectified. To guarantee that the portable system works perfectly, observe the following points.

## 7.1 Visual inspection

---

Look at the portable system closely. The portable system may not be damaged in any way

- which could result in parts coming loose
- which may restrict its function or
- which may cause injuries.

If you find any damage relating to one of these points, the portable system must be repaired or replaced.

## 7.2 Check catch

---

1. Hang the portable system in the wall mounting and allow it to engage.
  - The portable system must be located firmly in the wall mounting.
  - The red control pin of the catch cross-bar must be flush with the top edge of the catch housing.
2. Check unlatching next. To do this, reach through the handle and pull the catch cross-bar upwards.
  - The catch cross-bar must be easy to move.
  - The portable system may not catch on removal.

## 7.3 Check protective bag (if present)

---

### **Check push-buttons**

Open and close the push-buttons.

- The push-buttons may not pop open again.

### **Check Velcro tabs**

Open the Velcro tabs.

- Both types of tape should be free of contamination.

### **Check magnetic closures**

Open and close the magnetic closures.

- The magnetic closures may not pop open again.

## **8. Servicing**

---

The portable system and the wall mounting require no servicing, but observe the intervals for regular function checks on the attached parts. Look these up in the associated instructions for use if appropriate.

We recommend having maintenance measures such as servicing and repair work carried out by the manufacturer, WEINMANN Emergency, or by specialist staff.

### **8.1 Disposal**

---

#### **LIFE-BASE 1 NG without charger interface**

#### **BASE-STATION 1 NG without charger interface**

#### **BASE-STATION Mini II without charger interface**

Do not dispose of the device in domestic waste. To dispose of the device properly, contact a licensed, certified disposal company. You can find an address from your environment officer or from your municipal authorities. The device packaging (cardboard and inserts) can be disposed of in paper recycling facilities.

#### **LIFE-BASE 1 NG/LIFE-BASE 1 NG XL/LIFE-BASE 1 NG XS with charger interface**

#### **BASE-STATION 1 NG with charger interface**

#### **BASE-STATION Mini II with charger interface**

Do not dispose of the device in domestic waste. To dispose of the device properly, contact a licensed, certified electronic scrap disposal merchant. You can find an address from your environment officer or from your municipal authorities. The device packaging (cardboard and inserts) can be disposed of in paper recycling facilities.

# 9. Scope of supply

---

## 9.1 General information

---

The following applies to the portable system

- The portable system is supplied fully assembled depending on the equipment variant.
- Each portable system is supplied with these instructions for use.

The BASE-STATION 1 NG and BASE-STATION Mini II wall mountings are not included in the scope of supply and should be ordered separately if required (see “9.3 Accessories and spare parts” on page 29).

## 9.2 Standard scope of supply

---

Description	Order no.
<b>LIFE-BASE 1 NG without charger interface and</b>	
MEDUMAT Transport without CO <sub>2</sub> measurement	WM 9630
MEDUMAT Transport with CO <sub>2</sub> measurement	WM 9635
<b>LIFE-BASE 1 NG with charger interface and</b>	
MEDUMAT Transport without CO <sub>2</sub> measurement	WM 9620
MEDUMAT Transport with CO <sub>2</sub> measurement	WM 9625
<b>LIFE-BASE 1 NG XL with charger interface and</b>	
MEDUCORE Standard Advanced	WM 9720
MEDUCORE Standard Advanced plus MEDUMAT Easy CPR	WM 9725
MEDUCORE Standard Advanced plus MEDUMAT Standard a	WM 9730
MEDUCORE Standard Pro	WM 9755
MEDUCORE Standard Pro plus MEDUMAT Easy CPR	WM 9760
MEDUCORE Standard Pro plus MEDUMAT Standard a	WM 9765
MEDUCORE Standard Basic	WM 9735

<b>Description</b>	<b>Order no.</b>
MEDUCORE Standard Basic plus MEDUMAT Easy CPR	WM 9740
MEDUCORE Standard Basic plus MEDUMAT Standard a	WM 9745
MEDUMAT Standard <sup>2</sup>	WM 9400
MEDUMAT Standard <sup>2</sup> plus MEDUCORE Standard	WM 9405
<b>LIFE-BASE 1 NG XS with charger interface and</b>	
MEDUMAT Standard <sup>2</sup>	WM 9870
MEDUCORE Standard Basic	WM 9882
MEDUCORE Standard Advanced	WM 9883
MEDUCORE Standard Pro	WM 9880

## 9.3 Accessories and spare parts

---

You can order accessories and spare parts separately if required. A current list of accessories and spare parts can be ordered on the Internet at [www.weinmann-emergency.de](http://www.weinmann-emergency.de) or via your dealer.

# 10. Technical data

---

	<b>LIFE-BASE 1 NG</b>
Product class to 93/42/EEC	Class I
Standards applied	EN 1789
Dimensions (without protective bag) W x H x D in mm	415 x 240 x 140 (LIFE-BASE 1 NG) 427 x 240 x 140 (LIFE-BASE 1 NG XL) 220 x 230 x 136 (LIFE-BASE 1 NG XS)
Weight	Empty: LIFE-BASE 1 NG: approx. 2.0 kg LIFE-BASE 1 NG XL: approx. 2.2 kg LIFE-BASE 1 NG XS: approx. 1.0 kg With equipment: LIFE-BASE 1 NG and LIFE-BASE 1 NG XL: max. 7.0 kg LIFE-BASE 1 NG XS: max. 5.0 kg
Maximum load for carrying strap	5 kg
Temperature range – operation – storage	–18 °C to +60 °C –40 °C to +70 °C
<b>Note</b> The temperature range may be limited by any devices fitted!	
Electrical rating for charger interface	12 V DC, $I_{\max} = 8 \text{ A}$

 The right to make design modifications is reserved

# 11. Warranty

---

WEINMANN Emergency gives the customer a limited manufacturer warranty on new genuine WEINMANN Emergency products and any replacement part fitted by WEINMANN Emergency in accordance with the warranty conditions applicable to the product in question and in accordance with the warranty periods from date of purchase as listed below. The warranty conditions can be downloaded from [www.weinmann-emergency.de](http://www.weinmann-emergency.de) on the Internet. We can also send you the warranty conditions on request.

In the event of a claim under warranty, contact your specialist dealer.

Product	Warranty period
WEINMANN Emergency devices including accessories (exception: masks) for oxygen and emergency medicine	2 years
Masks including accessories, rechargeable battery, batteries (unless there is information to the contrary in the technical documentation), sensors, tube systems	6 months
Disposable products	None

## 12. Declaration of conformity

---

WEINMANN Emergency Medical Technology GmbH + Co. KG declares herewith that the product complies fully with the respective regulations of the Medical Device Directive 93/42/EEC. The unabridged text of the Declaration of Conformity can be found on our website at [www.weinmann-emergency.de](http://www.weinmann-emergency.de)

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